

SIEMENS

POWERMOBIL

SP

System Manual

Installation Instructions

Unpacking the Monitor Trolley

© Siemens AG 1999

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Register 3

Print No.: SPR2-240.812.02.02.02

Replaces: SPR2-240.812.02.01.02

English

Doc. Gen. Date: 02.01

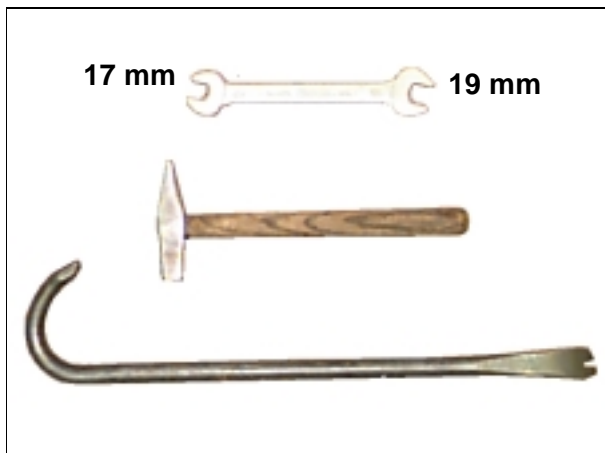


Fig. 1 Required tools



Fig. 2 ① Ramp, Pno. 55 63 783 G5464

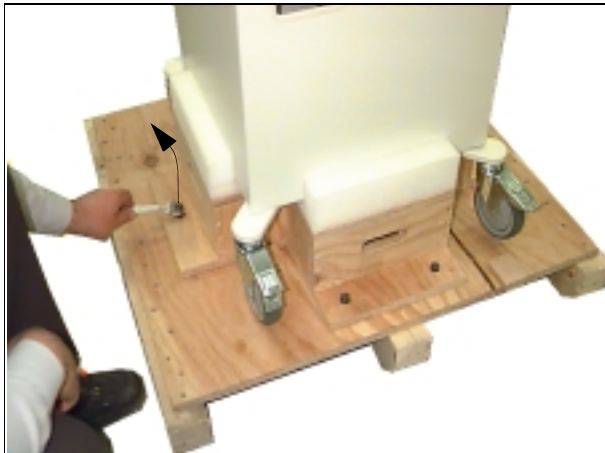


Fig. 3 Remove screws from side support blocks



Fig. 4 Remove first support block (right)



Fig. 5 Remove second support block (left)

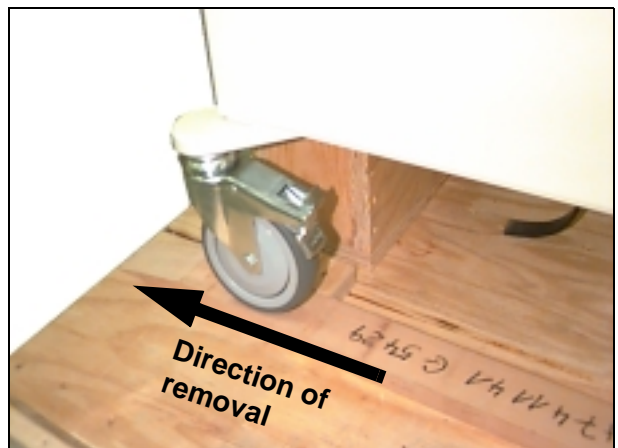


Fig. 6 Line up castors / Lock brakes

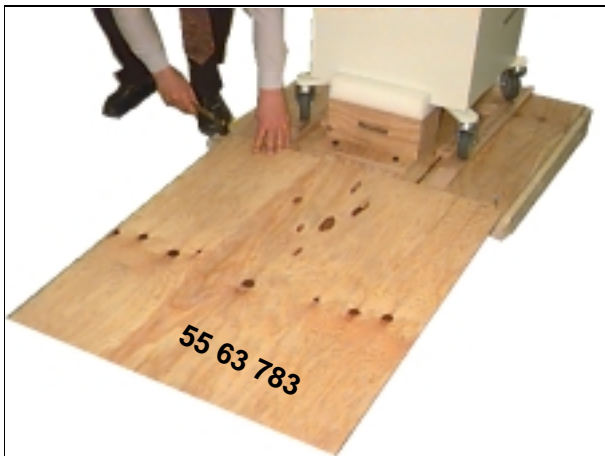


Fig. 7 Fasten ramp 55 63 783



Fig. 8 Remove screws from front and rear of support blocks



Fig. 9 Remove front and rear support blocks



Fig. 10 Unlock all castor brakes / Roll unit off of pallet

Concluding work

CAUTION

In order to avoid damages caused by relative humidity, do not put the POWERMOBIL into operation immediately. In extreme cases, it can take up to 24 hours before the POWERMOBIL should be turned on. (transport temperature -20° Celsius - very high humidity in the customer's clinical rooms).

- The basic unit should be unpacked if this has not been done.
- Prior to start-up, let the system adjust to its environment.
- Perform the system start-up according to the instructions outlined in the POWERMOBIL Start-up manual.

TD PS 24 Arnold
A&D PT 6 Dr. McMinn

Changes to previous version

Page 1	Cover page, revision level changed
Page 3	Fig. 7: Number changed